

Motion Sensor Setup & Install

Aloe Care



Help Is Here

Reach out to us for assistance with setup or answers to general questions; one of our team members will be happy to help.

Email, or speak to us directly—

support@aloecare.com


1-888-ALOECARE (1-888-256-3227)

support.aloecare.com

Contents

1. Your Motion Sensors	04
2. Insert the Batteries	06
3. Close the Casing	08
4. Pairing the Device	10
5. Tabletop Placement (best)	12
6. Wall Mounting	14

One Set

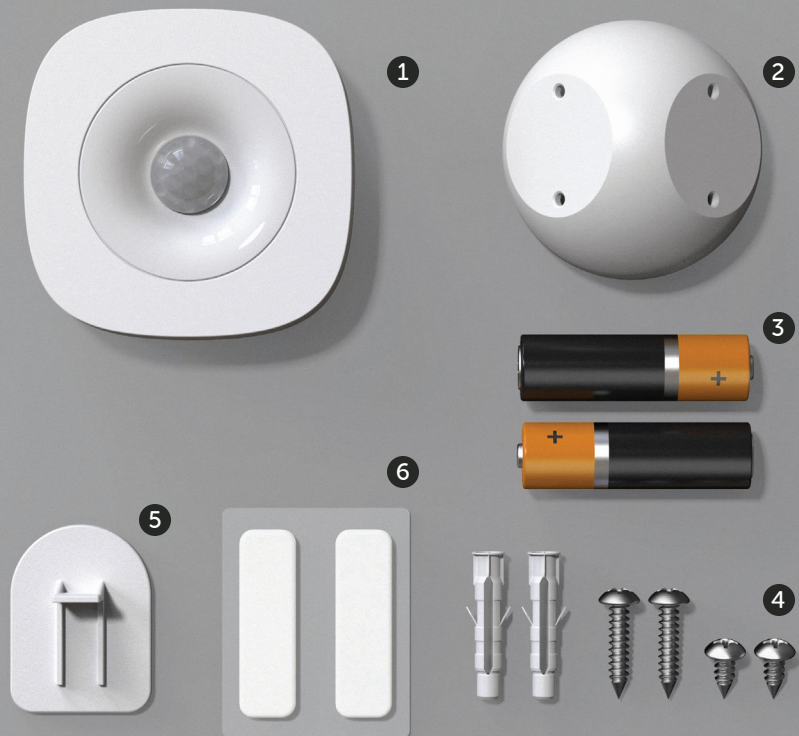


1

Meet Your Motion Sensors

Unpack the two sets of Motion Sensors and lay the items out on a table for easy setup.

Check your items—

1. Motion Sensor
2. Corner Mount
3. AA Batteries
4. Corner Mount Accessories
5. Mounting Tape
6. Device Stand

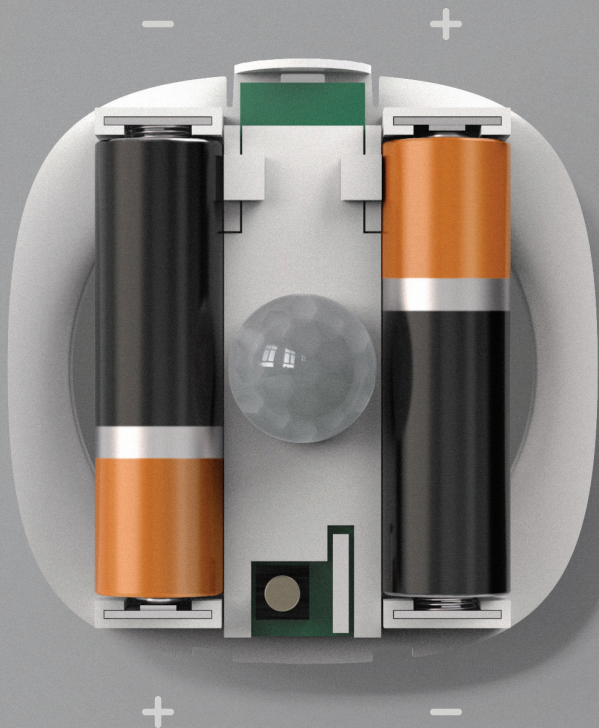

2

Insert the Batteries

To start, press the circular tab on the outside edge of the Motion Sensor casing and separate the front and back panel. Then, insert the two AA batteries into the back panel, making sure the negative poles touch the springs.

TIP Once the batteries are inserted properly, the sensor will begin flashing red.




3

Close the Casing

To close the Motion Sensor, line up the teeth at the bottom of both panels (they fit together), then press the casing closed.

TIP The two panels will “click” together when assembled correctly.

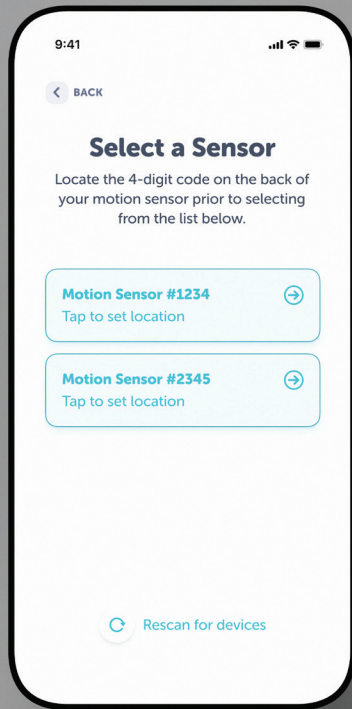

4

Pairing the Device

To pair the Motion Sensor to the Smart Hub, open the Aloe Care app on your mobile device and navigate to the Device Management menu.

Referring to the four-character code on the back of your Motion Sensor, select the corresponding card from the app. Follow the on-screen instructions to finish pairing your device. Repeat process for additional sensors.

TIP Make sure you're sensors are near the Smart Hub when pairing.



5

Tabletop Placement (recommended)

For tabletop placement, insert the stand into the slot on the bottom of the device. Then, place the sensors in commonly used rooms not already occupied by a smart sensor, and where they will have a clear view of the room.

Do not—

- ✗ Block motion sensors.
- ✗ Place near a water fixture.
- ✗ Place in direct sunlight.
- ✗ Place near heating, or cooling.



For installation instructions using the corner mount, go to our online support page at support.aloecare.com.

6

Mount It on a Wall

Find a spot on a wall, where the Motion Sensor will have a clear view of the room. Apply both pieces of the double-sided tape to the back panel, then firmly stick the device on the spot you've chosen.

For proper function—

- ✗ Do not mount on porous or textured walls, the tape will not stick well.
- ✗ Place at least 14 feet from a fireplace.
- ✗ Do not mount on a ceiling using the tape.

Aloe Care Health, Inc.
Copyright © 2020 Aloe Care Health