


Care Button Setup

Aloe Care

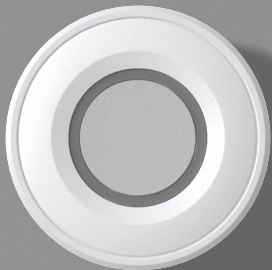

Meet Your New Care Button

The Care Button is a waterproof accessory device that extends the Smart Hub's reach to other areas of the home and property. When pressed, the Care Button pings the Smart Hub and sends a notification to the care network.

This device operates as an alert beacon and does not contain a speaker or a microphone.



1

Activate It

An Aloe Care team member must activate the Care Button. Call the number below to get started.

TIP Have the Care Button box handy for the call—you will need to reference a number on it for activation.

1-888-ALOE CARE
(1-888-256-3772)



2

Wear It

Snap it into the lanyard to wear it like a necklace, or into the wristband to wear it like a watch.

3

Press to Notify

Press the button to get assistance from your care network. When pressed, the button's ring light turns green.

Aloe Care Health, Inc.
Copyright © 2020 Aloe Care Health